



HOW TO ADD RACE ALIGNERS TO YOUR 3SHAPE SCANNER




- 1 Go to 3Shape Communicate portal <https://portal.3shapecommunicate.com/login>
- 2 Click 'Sign in' button

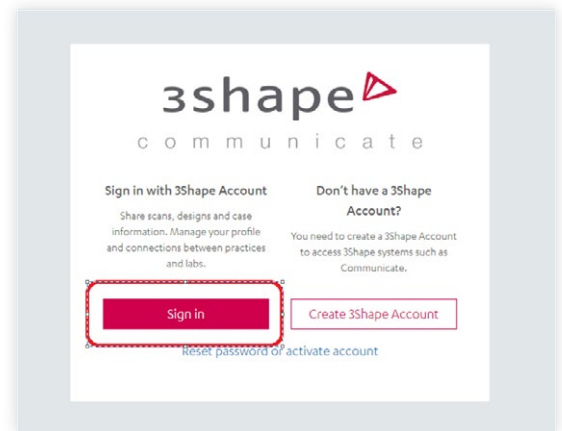


- 3 Please use your dentist's communicate account to login.

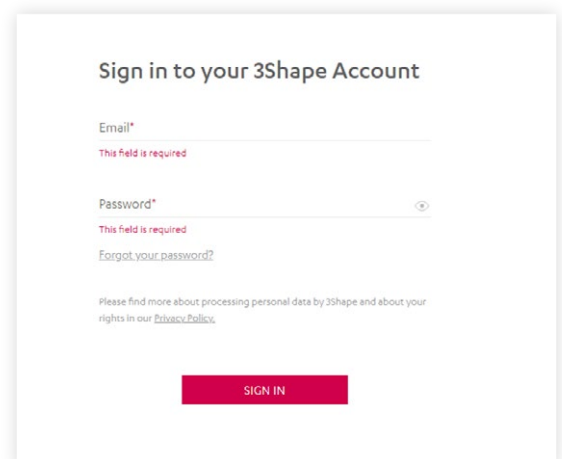


- 4 Once you logged in the 3Shape Communicate portal, please choose the tab of **'Connections'** and click **'Add Connection'** button.

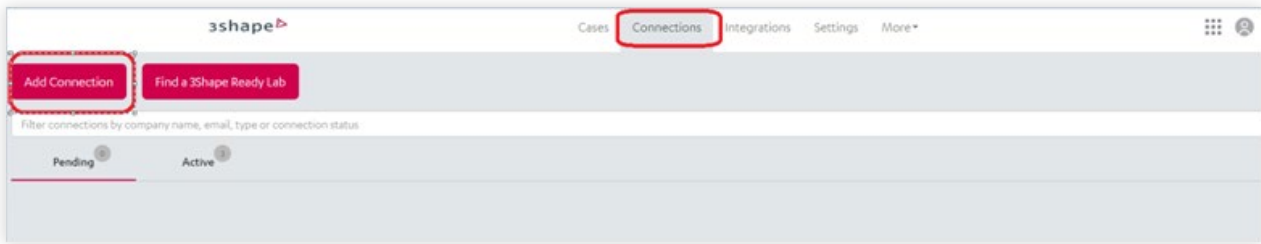

- 5 With the popped up search bar, please input Race Aligners account **'racealigners@racedental.com.au'**.

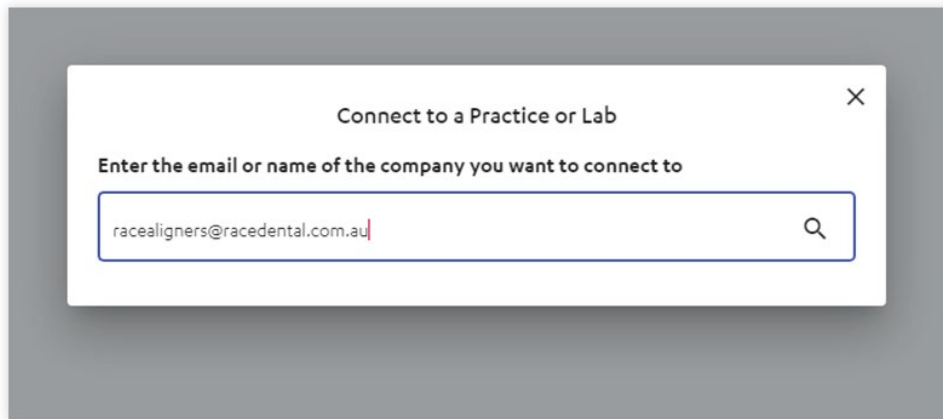


- 6 Then you will be able to see our Aligners account and please click **'Connect'** button.

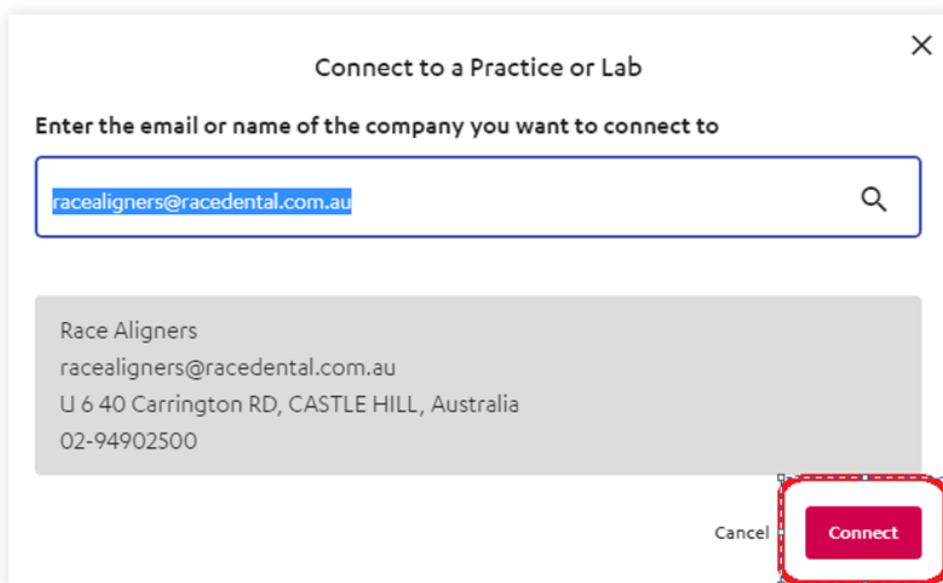


HOW TO SEND YOUR **ALIGNER** **CASE** FROM YOUR **3SHAPE SCANNER**




- 1 Set up a new patient case. Select **Race Aligners**.

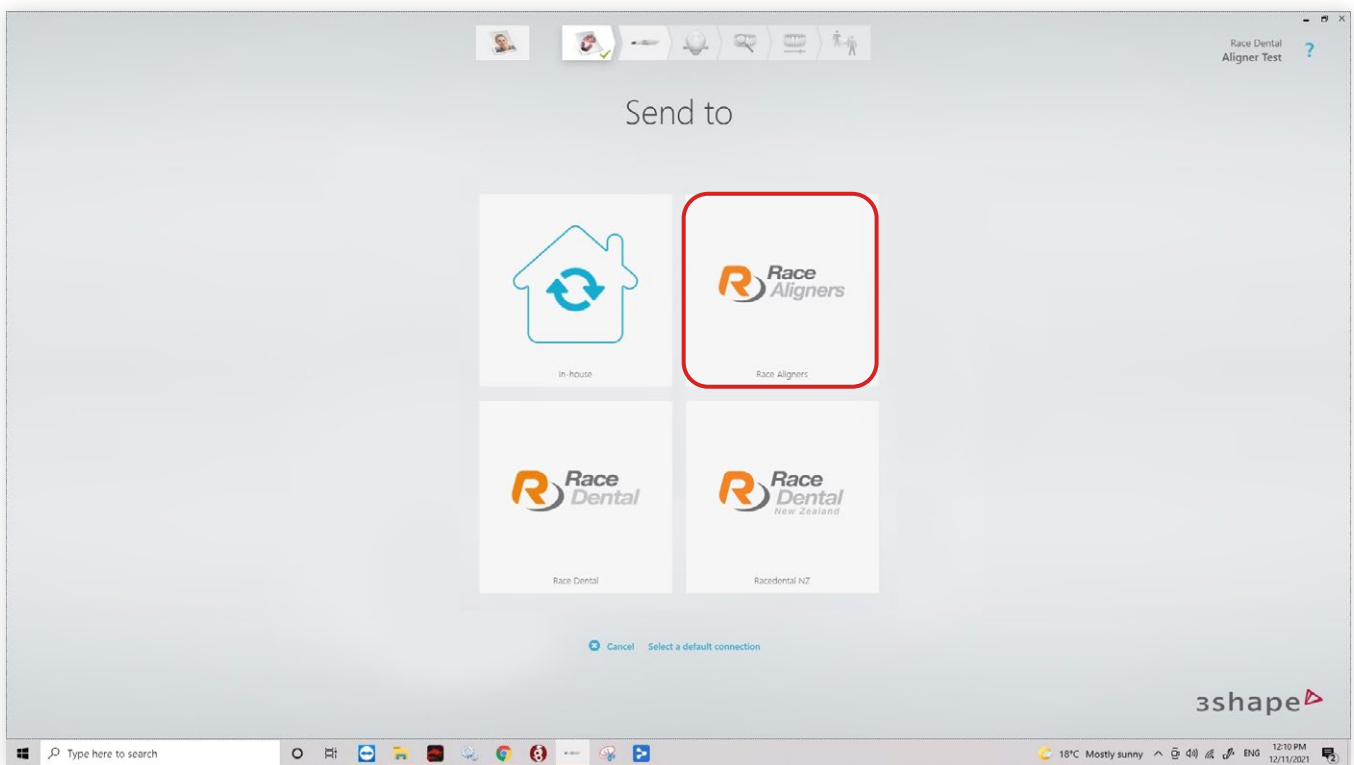

2 In the prescription select **Orthodontics**.

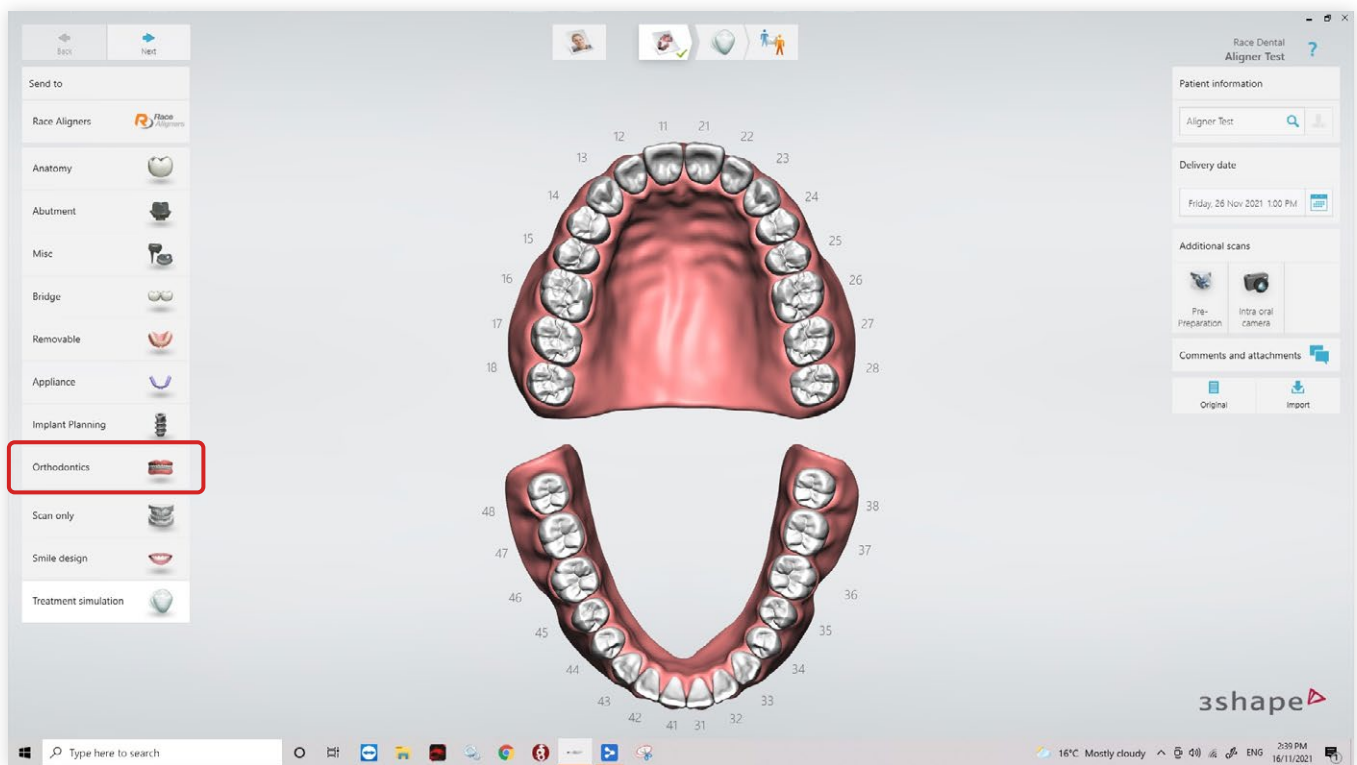


3 Then click **clear aligners**.

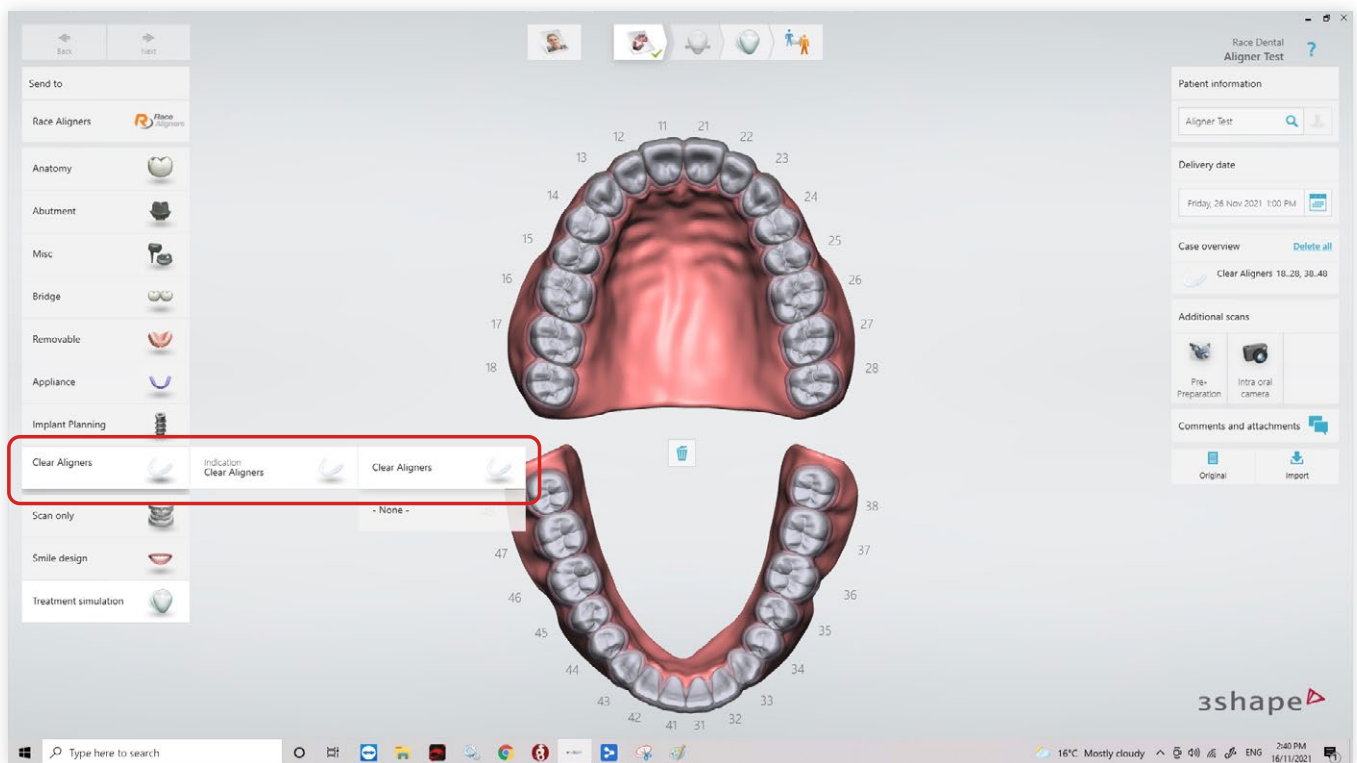


You may include any comments or attachments on the right.

4 Make sure to set a delivery date.

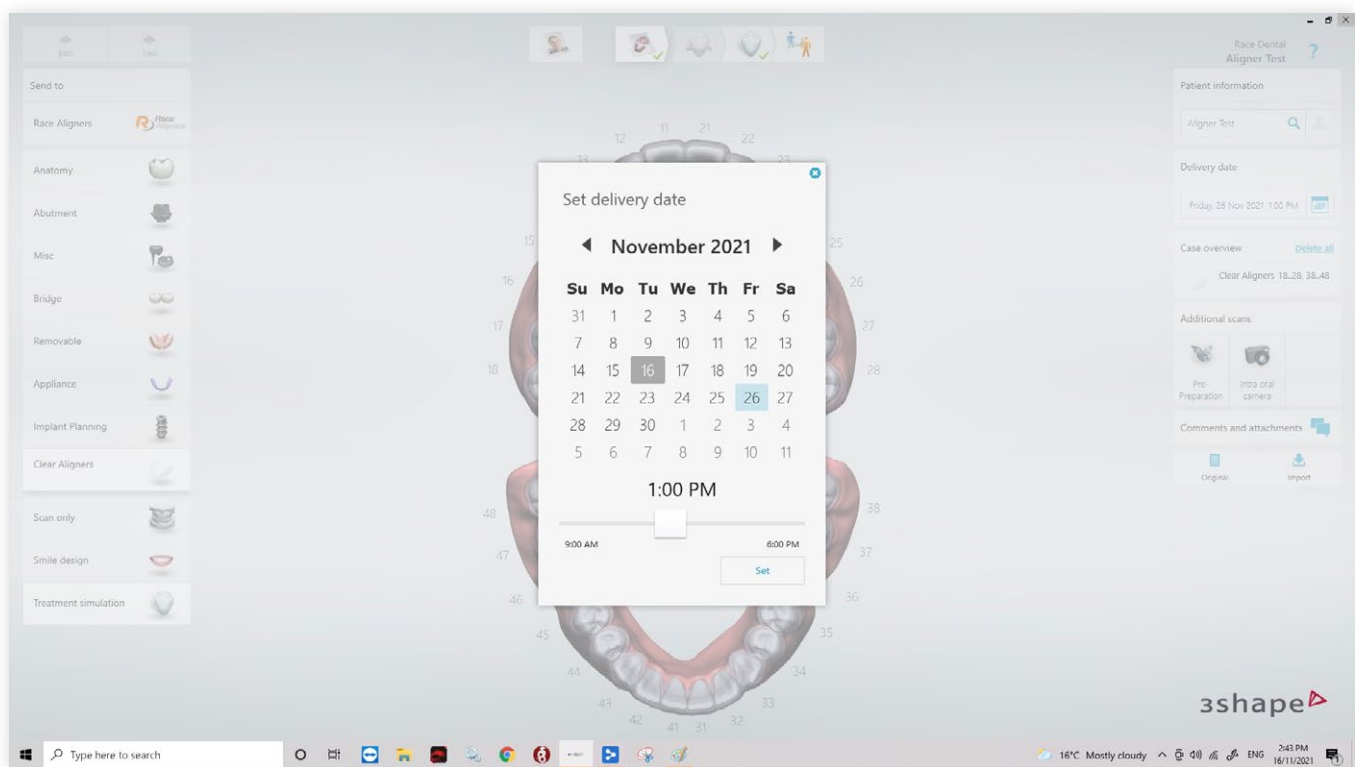

5 Once full mouth scan is completed and processed. Send Case.

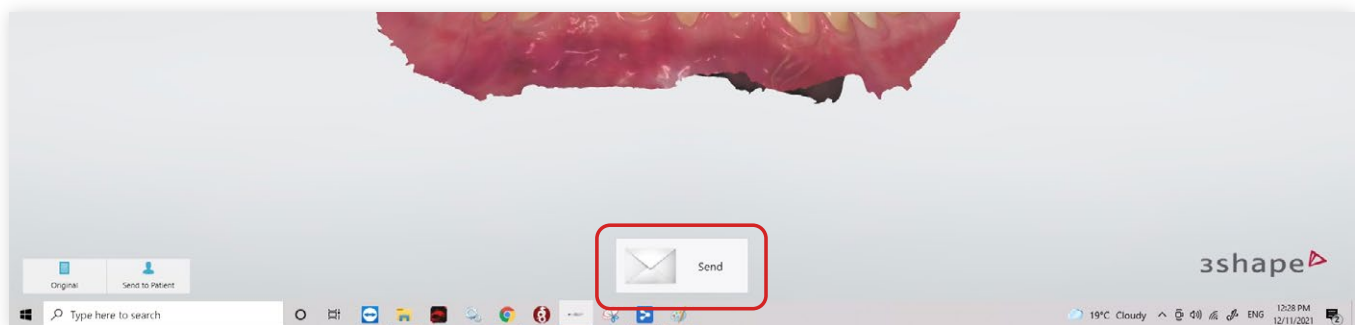

6 Complete the Race Aligners Rx Form and send to customersupport@racedental.co.nz

Race Aligners RX Form



Doctor Information

First Name Last Name

Shipping Address

Patient Details

First Name Last Name

Chart Number Sex Male Female

Case Details

Shade

Scanner Traditional impression 3M 3Shape Carestream
 iTero Medit Other Digital

Submission Date

Requested Return Date

Arch Upper Lower Upper & lower

Patient's Mobile
**required for dental monitoring*

Patient's Date of Birth
**required for dental monitoring*

Patient's Email
**required for dental monitoring*

Allow IPR Yes No

Allow Attachment Yes No

Race Aligners RX Form



Indicate Extractions

**if applicable*

18	17	16	15	14	13	12	11	21	22	23	24	25	26	27	28
<hr/>								<hr/>							
48	47	46	45	44	43	42	41	31	32	33	34	35	36	37	38

Ankylosis/Implant

**tooth that cannot be moved
- if applicable*

18	17	16	15	14	13	12	11	21	22	23	24	25	26	27	28
<hr/>								<hr/>							
48	47	46	45	44	43	42	41	31	32	33	34	35	36	37	38

AP Relation - Left

Maintain Improve Canine Relationship

AP Relation - Right

Maintain Improve Canine Relationship

Overjet

Maintain Improve

Overbite

Maintain Improve

Comment/Further Specification